


## BULK LIQUID SUPPLY

- Continuous loading
- Increased walk-away time
- Prevention of system interruptions
- User-friendly container exchange
- Prevention of spillage

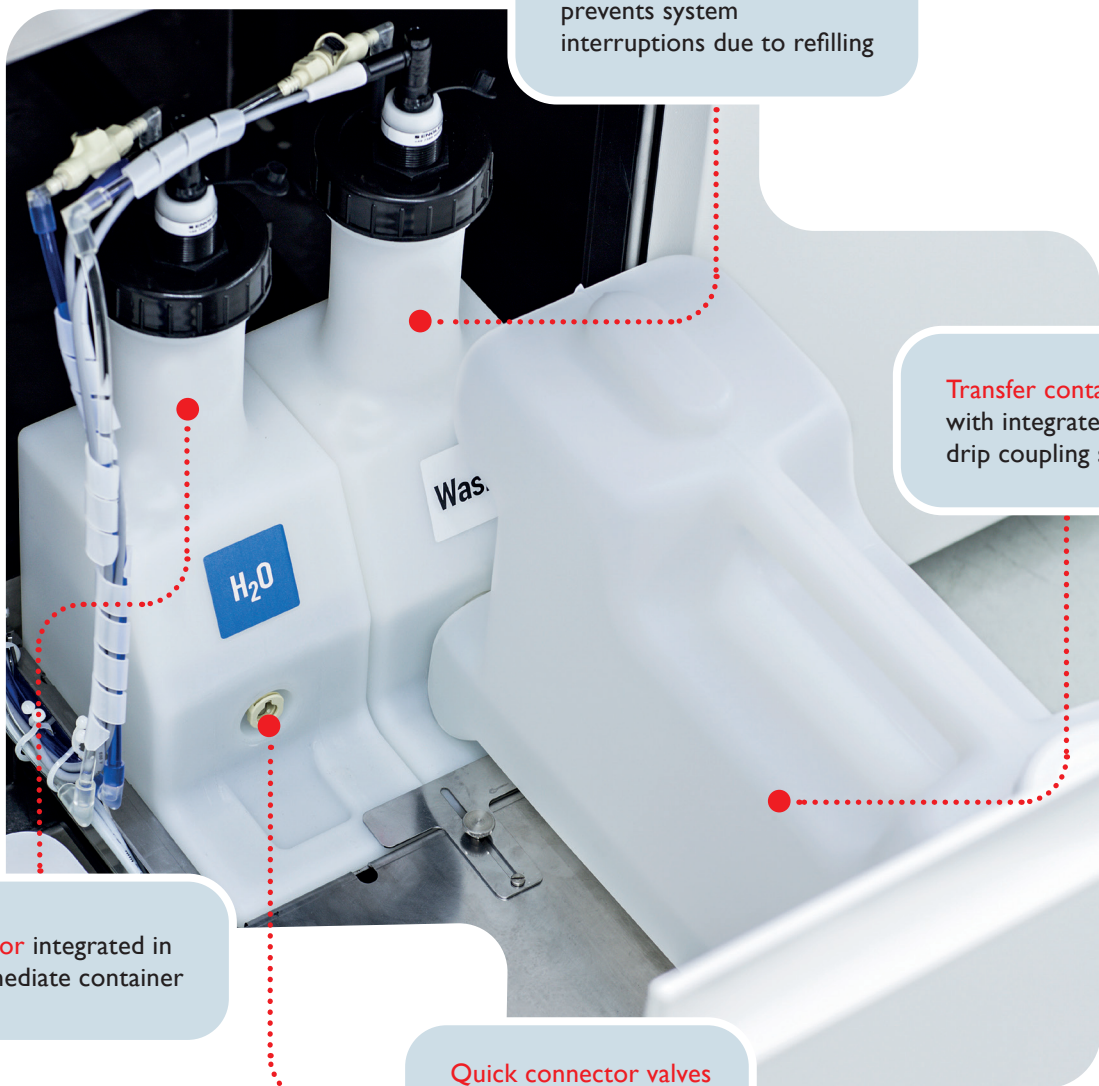

Intermediate container prevents system interruptions due to refilling

Transfer container with integrated anti-drip coupling system

Level sensor integrated in the intermediate container

Quick connector valves enable easy container exchange

Intermediate containers with level detection and quick connector systems enable refilling of system liquids without interruptions and spillage. There is no need to switch tubes or level detection from one container to another.

In the transfer containers integrated anti-drip quick connectors prevent liquid spillage.

Bulk Liquid Supply is perfectly suited for liquids required in high volumes, e.g. system liquids like DI water and wash buffer.

## SPECIFICATIONS

<b>Volume transfer container</b>	Up to 10 l
<b>Volume intermediate container</b>	Up to 3.8 l
<b>Weight (empty module)</b>	Approx. 10 kg
<b>Refill opening size</b>	80 mm
<b>Dimensions [H W D]</b>	393 368 749 mm
<b>Wetted components</b>	PP: CPC connectors, level sensor; HDPE: container

The herein provided information and specifications by STRATEC SE are provided on an "as-is" basis and may be subject to changes, omissions and errors. STRATEC SE shall have no liability for the accuracy and completeness of these information and specifications and cannot be held liable for any representations or warranties, expressed or implied, and any third party claims or losses of any damages thereof. STRATEC SE reserves the right to make any changes to the information and specifications at any time and without notice. The applications and products herein may not be available in all markets.

STRATEC is a registered trademark. © 2022, STRATEC SE, Germany, all rights reserved.